

Using attentions of heating pads and adapters

Please reading the attentions detailedly as follows

- Please ensure the normal operation of the product according to the special transformer and power outlet made by the manufacturer.
- Please don't put the heating pads off the sunshine directly or near the inflammable goods.
- Please don't disassemble by yourself so as not to be electrocuted.
- Please don't scrub with acidic or alkaline detergent.
- Please don't put product on high temperature or high humidity place to use.
- Please disconnect the power supply and unplug the power adapter when the product is not using or person is leaving .
- The power line can not be knotted, bunched together etc. and do not use it when the husk is damaged.
- Please don't use the multi-function tranfer socket.
- When the following situations are encountered,please unplug power supply line and stop use it.

*The heating pads and power adapters are affected with damp or rain wet.


* The smoke, sparks and other safety factors occur on the line.

*The destory or damage of the power line and power socket.

(If there are have any problems please contact the after-sale customer service staff.)

- Pregnant women with caution
- We are not responsible for the failure of the product caused by the disassembly of the product by yourselves.

Note: In order to your safety and avoid accidents, please follow the above attentions to use. We are not responsible for any damage or accident caused by using heating pads and power adapters without following the instructions above.

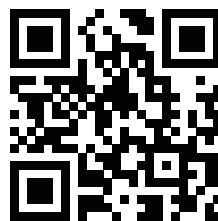

Phototherapy Pad

Phototherapy Pad

Instruction manual



Please read the instructions carefully before you use the product

Shenzhen Guangyang Zhongkang Technology Co.,Ltd

+86-13631598551

www.suyzeko.com

3rd Floor, Building N, Innovation Silicon Valley Industrial Park,
17 Jingtian Road, Guanhu Street, Longhua District, Shenzhen

*Please read this instruction carefully then use it correctly and keep it properly

Packing list



LED photon pad	1	Power line	1
Velcro	1	Manual	1
USB charging head	1	Pabric packing bag	1

Product function


- * LED photon pad is a kind of health care product that promotes deep subcutaneous heat through deep infrared heat, thereby promoting blood circulation, relieving soreness, eliminating fatigue and preventing disease. Suitable for rheumatism, rheumatoid arthritis, bone hyperplasia, frozen shoulder, cervical spondylosis, low back pain, stomach cold, hand and foot paralysis.
- * Red light-Trigger microcirculation increase.
- * Infrared-Reach deeper tissues, such as muscles, tendons, and ligaments.

Operation Instructions

①Take the photon pad and power cable out from the inner of package box separately and then connect the pad line to the power cable with the side of control box.

②Connect the line with power supply via the USB port, powerbank and power adapter are available.

③Finally the installation is complete. Turn on the button to make it work.

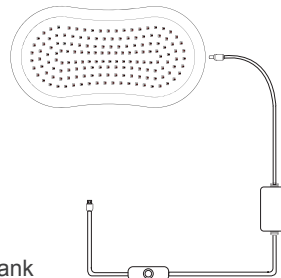

* Note:

After connecting to the power supply, it will light up when plugged in, and will stop automatically after 30 minutes of continuous use.
The small black box has a cycle function with frequency combination.

Basic parameter

Product name	Phototherapy Pad	Infrared wavelength	850nm
LED quantity	152pcs	Red light wavelength	660nm
Total power	6W	Output power	5V --- 2A
Input voltage	100-240V	Input current	2A

After-sale service guarantee card

- Since the purchase of this product, under the normal use of the guarantee of one year free service. If beyond the guarantee period will charge a certain cost of parts.
- If you need the replace parts service, please according to our service spare parts maintenance price list to charge fee.
- When the product returns to the factory for repair, it should be properly packed when it is sent to repair, so as not to result damage in the course of delivery.
- Under following conditions, although in the gurantee period, it can't repair without fee.
 - 1.A breakdown or damage caused by natural disasters, rodent disasters, and improper use.
 - 2.A breakdown or damage caused by delivery, movement, and improper preservation after purchase.
 - 3.The information record is incomplete, does not fill in the purchase date and purchase shop stamp.
 - 4.Special environment is not guarantee, such as wet, sulfur areas, coastal high salt areas and ships or other corrosion-prone areas.
 - 5.The user who didn't follow the manual.